


ageLOC®
WellSpa iO™

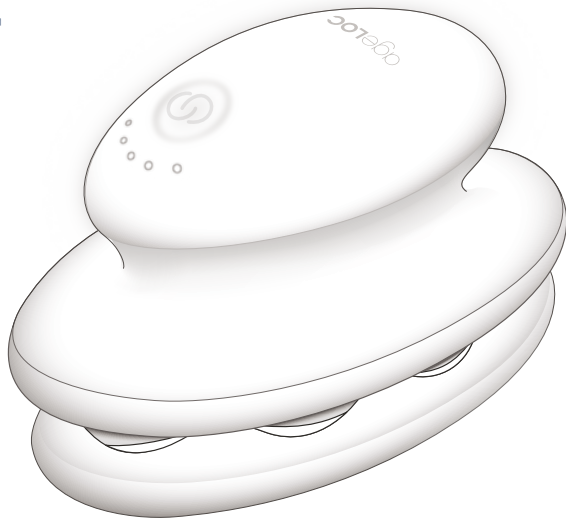

USER'S MANUAL



TABLE OF CONTENTS

System Contents and Components	1–2	Maintenance	13
Safety Instructions, Warnings, and Precautions	3	Replacement and Warranty Information	13
Charging Your ageLOC® WellSpa iO™	6	Patents	13
Connecting Your ageLOC WellSpa iO to the Nu Skin Vera® App	7	Technical and Regulatory Information	14
Pairing Your ageLOC WellSpa iO	7	Definition of Symbols	14
Using Your ageLOC WellSpa iO	8		
Cleaning and Caring for Your ageLOC WellSpa iO	10		
Troubleshooting	11		

SYSTEM CONTENTS


ageLOC® WellSpa iO™ device

Quick Start Guide

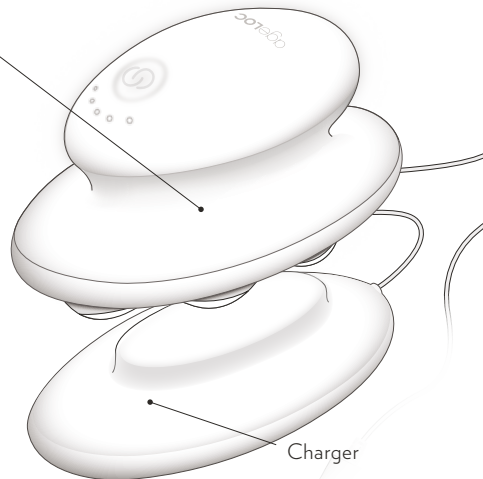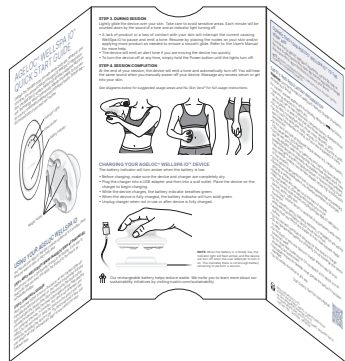

USB charging adapter not included.


SYSTEM COMPONENTS

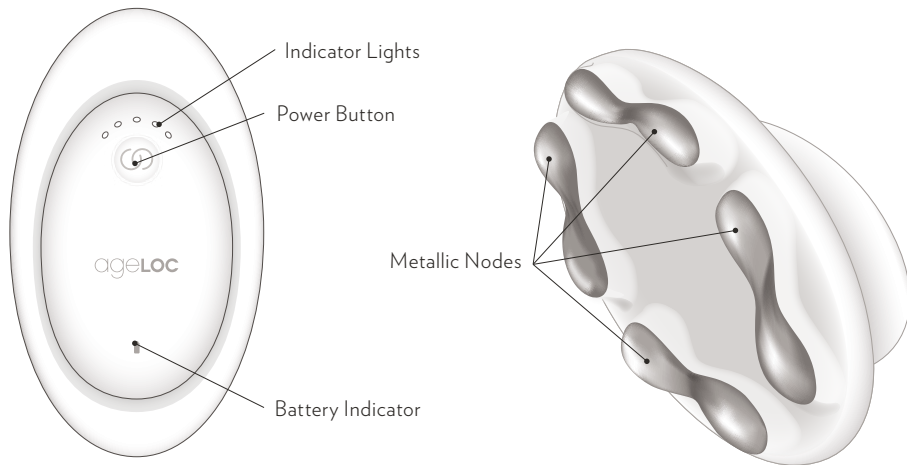

SAFETY INSTRUCTIONS, WARNINGS, AND PRECAUTIONS

ageLOC® WellSpa iO™ and the ageLOC WellSpa iO Wireless Charger are tested at pressures matching 4000m in altitude. This appliance contains batteries that are non-replaceable.

WARNINGS

 **WARNING** Only use the ageLOC WellSpa iO Wireless Charger provided to recharge your device.

- The ageLOC WellSpa iO device is only for use on healthy skin.
- Do not apply to head/face, neck, throat, chest, décolleté, or groin. Do not use over areas of skin that lack normal sensation as it could result in irritation.
- Do not use if you have open sores, wounds, or rashes, are prone to pityriasis rosea, have allergies to metal, or if you have problematic or overly sensitive skin.
- Do not use over swollen, red, infected, or inflamed areas; or over skin eruptions (e.g., phlebitis, thrombophlebitis, varicose veins).
- Do not use over, or in proximity to, cancerous lesions or diseased tissues.
- Do not touch another person, or electrically live parts/items while using the device
- Do not use for longer periods of time or more frequently than recommended.
- Do not operate while bathing or showering. Do not submerge your device in water.
- Do not use if you have an implanted cardiac device (e.g., cardiac monitors, ICDs, pacemakers, ECG alarms).
- Consult with your physician or healthcare provider if you have sensitive skin, are epileptic, wear body piercings, have a metal implant, or are ill. If you experience irritation or prolonged stinging, redness, discontinue use immediately and seek medical advice, if necessary.
- Consult with your physician or healthcare provider before using the device if you have bleeding tendencies, bruise easily, are pregnant or may become pregnant, have undergone a recent surgical procedure.

SAFETY INSTRUCTIONS, WARNINGS, AND PRECAUTIONS (Cont.)

- Consult with your physician or healthcare provider before using the device if you have any ongoing or previous medical conditions, such as heart disease, stroke, high blood pressure or liver disorders.
- Keep the device and charger out of the reach of children.
- Remove jewelry or other metal adornment from the surface of the skin prior to use.
- Do not use with other devices, including multiple, ageLOC® WellSpa iO™ devices, at the same time.
- ageLOC WellSpa iO is for individual use only.
- Charge only with the ageLOC WellSpa iO Wireless Charger. No other charging accessories should be used with the device.
- Always place your ageLOC WellSpa iO device and geLOC WellSpa iO Wireless Charger on a dry, heat-resistant, stable, flat surface, especially when charging.
- Do not keep the charger plugged in when not in use or after device has reached full charge.
- Do not use harsh chemicals or abrasives on your ageLOC WellSpa iO device or ageLOC WellSpa iO Wireless Charger.
- To reduce the risk of shock, irritation, fire, or injury, do not place or store your device or charger where it may fall or be pulled into water.
- Do not place or store the device where it may fall or be pulled into a bathtub, sink, or toilet.
- Do not store near a heat source, such as a radiator, fire, or heat vent. Do not leave in a hot vehicle.
- Do not leave the ageLOC WellSpa iO device in an extremely cold environment for prolonged periods of time.
- Do not use the device to apply electrical stimulation while driving, operating machinery, or engaging in any activity in which electrical stimulation can put you at risk of injury.
- Periodically inspect for damage; never operate the device if it has been damaged. In the case of any device defect or if the device gets damaged, contact Nu Skin Support Services.
- Do not modify the device in any way.

TEMPERATURE

Ideal device storage and charging is to occur at room temperature, between 20°C and 24°C (68°F to 75°F).

Only store the device in environments between 5°C and 45°C (41°F to 113°F).

Storage and charging outside of these ranges will adversely affect the battery, metallic nodes, and operation of the device.


Very high temperatures or hot environments, (e.g., above 50°C/122°F, direct sunlight, within vehicles in extremely hot conditions, etc.) can result in overheating and other catastrophic conditions that severely impact the performance and life of the device.

Very low temperatures or freezing environments can result in permanent damage to the metallic nodes and operation.

CHARGING YOUR AGELOC® WELLSPA iO™ DEVICE

The battery indicator will turn amber when the battery is low.

- Before charging, make sure the device and charger are completely dry.
- Plug the charger into a USB adapter and then into a wall outlet. Place the device on the charger to begin charging.
- While the device charges, the battery indicator breathes green.
- When the device is fully charged, the battery indicator will turn solid green.
- Do not keep the charger plugged in when not in use or after device has reached full charge.

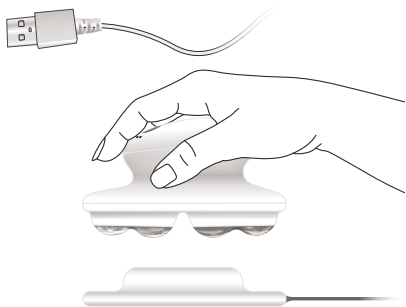

NOTE

When the battery is critically low, the battery indicator will flash amber, and the device will turn off when the user attempts to turn it on. This indicates there is not enough battery remaining to perform a session.

If there are no indicator lights, charge the device. Indicator lights will turn on once the device is charged.



Our rechargeable battery helps reduce waste. We invite you to learn more about our sustainability initiatives by visiting nuskin.com/sustainability.

CONNECTING AGELOC® WELLSPA iO™ TO THE NU SKIN VERA® APP

To unlock the full ageLOC WellSpa iO experience, scan the QR code to download the Nu Skin Vera® app from the App Store or Google Play. Follow the app instructions to connect to your device.



PAIRING YOUR AGELOC® WELLSPA iO™

1. Download the Nu Skin Vera® app.
2. Ensure your phone's Bluetooth®* is enabled.
3. Sign in to Nu Skin Vera.
 - Already have a Nu Skin account? Use it to sign in.
 - New to Nu Skin? Sign up—it's easy!
4. Tap "Devices" in the navigation bar.
5. Tap "Pair New Device"
6. If the device is off, press the Power button to turn the device on.
7. Tap "Connect."
8. The device will pair and you're ready to go.

* The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by NSE Products, Inc. is under license. Other trademarks and trade names are those of their respective owners.

App Store and the Apple logo are trademarks of Apple Inc., registered in the U.S. and other countries. Google Play is a trademark of Google LLC. For information about how Nu Skin collects, uses, and shares Personal Data please visit <https://www.nuskin.com/content/global-privacy.html>

USING YOUR AGELOC® WELLSPA iO™ DEVICE SYSTEM

You may use your ageLOC WellSpa iO device system up to 5 minutes, three times per week to daily, on each body area.

STEP 1. APPLY AGELOC BODY SERUM OR AGELOC BODY ACTIVATING GEL

Smooth your serum or gel liberally on the area to be treated, avoiding face, neck, chest, and groin. If you plan to treat more than one area, only apply the product to the area to be immediately treated.

STEP 2. STARTING A SESSION

Press the power button to turn on the device. All five indicator lights will illuminate the first time you power on. The number of illuminated indicator lights indicates the session time. Each time you press the Power button, you select a new time. Press the Power button again until you have selected your preferred time from 1–5 minutes. Your device will remember your setting for the next session. Place the device on the area of your skin where you've applied the serum or gel. The device will emit a tone, and your session will start automatically. Begin a new session for each treatment area.

STEP 3. DURING THE SESSION


Lightly glide the device over your skin. Take care to avoid sensitive areas. Each minute will be counted down by the sound of a tone and an indicator light turning off.

- A lack of product or a loss of contact with your skin will interrupt the current causing WellSpa iO to pause and emit a tone. Resume by placing the metallic nodes on your skin and/or applying more product as needed to ensure a smooth glide.
- The device will emit an alert tone if you are moving the device too quickly.
- To turn the device off at any time, simply press and hold the Power button until the lights turn off.

STEP 4. SESSION COMPLETION

At the end of your session, the device will emit a tone and automatically turn off. You will hear the same sound when you manually power off your device. Massage any excess serum or gel into your skin.

See diagrams below for suggested usage areas and Nu Skin Vera® for full usage instructions.

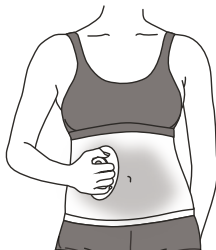

CLEANING AND CARING FOR YOUR AGELOC® WELLSPA iO™

Use slightly damp, lint-free wipe(s) to wipe away any serum, gel, or debris from the handheld device and charger after each use. Continue to wipe until visibly clean. Store in a cool, dry place. Do not immerse device in water or any liquid at any time. Using dry, non-linting wipe(s), thoroughly dry the external surfaces of the handheld device and charger.

TROUBLESHOOTING

LOSS OF CONDUCTIVITY DURING USE

If the metallic nodes lose contact with your skin, the current will be interrupted, causing WellSpa iO to pause, emit a tone, and the indicator lights to breathe.

PRO TIP: A lack of serum or gel may cause a loss of conductivity. Apply more product to your skin and ensure the device is in contact with both your skin and the product for optimal use.

MOVEMENT ALERTS DURING USE

The device will emit a sound if you are moving too quickly. The session will automatically resume when you resume the correct movement speed. Slowly glide the device over the targeted areas.

RESET BLUETOOTH® CONNECTION

If you have problems pairing with Nu Skin Vera® or if you need to unpair your device from your phone, conduct a Bluetooth reset.

1. Disconnect from Nu Skin Vera within the app by navigating to Device Settings > Forget this Device.

2. Disconnect your phone.

- **iPhone:** Navigate to your phone's Settings > Bluetooth > My Devices > NuSkinBDiO > Forget This Device.
- **Android:** Navigate to your phone's Settings > Connections > Bluetooth > Settings icon > Forget/Unpair.

3. Place WellSpa iO on the charger.

4. Press and hold the Power button for 5 seconds.

Wait for:

- Indicator lights 1-5 to flash once.

*The Bluetooth® word mark and logos are registered trademarks owned by Bluetooth SIG, Inc. and any use of such marks by NSE Products Inc. is under license. Other trademarks and trade names are those of their respective owners.

FACTORY RESET

To factory reset your ageLOC® WellSpa iO™:

1. Disconnect from Nu Skin Vera® within the app by navigating to Device Settings > Forget this Device.
2. Disconnect from your phone.
 - **iPhone:** Navigate to your phone's Settings > Bluetooth® > My Devices > NuSkinBDiO > Forget This Device.
 - **Android:** Navigate to your phone's Settings > Connections > Bluetooth > Settings icon > Forget/Unpair
3. Place WellSpa iO on the charger.
4. Press and hold the power button for 10 seconds.

Wait for:

- Indicator lights 1–5 to flash once.

Keep pressing and holding until:

- All indicator lights flash 3 times (at 10 seconds).

NOTE

ageLOC WellSpa iO automatically stores quality and usage information. When the device is factory reset, some device usage data will be retained for quality purposes. To view Nu Skin's privacy notice, visit:

-https://www.nuskin.com/en_US/corporate/privacy.html

DEVICE DISPOSAL

You must dispose of ageLOC® WellSpa iO properly according to local laws and regulations. The device contains electronic components and a lithium-ion battery, which requires disposal separate from household waste. When ageLOC WellSpa iO reaches its end of life, contact your local waste management department to learn about proper disposal and recycling options.



Dispose of batteries properly for your area.

MAINTENANCE

See the “Cleaning and Caring for Your ageLOC® WellSpa iO” section for routine maintenance.

Do not ship damaged devices by air.

SERVICING

Do not attempt to repair the device yourself, as this will void any warranty. There are no user-serviceable parts. Please see the warranty section for service.

REPLACEMENT AND WARRANTY INFORMATION

Limited One-Year Warranty: Nu Skin guarantees your device to be free from defects in materials and workmanship for a period of one year from the original date of consumer purchase. This warranty does not cover damage to the product due to misuse or accident, including dropping the device. If the product becomes defective within the one-year warranty period, please call your local Nu Skin Support Services to arrange for a replacement.

If your device has been discontinued, Nu Skin may replace your device with a new product which meets or exceeds the specifications and performance of your device. In such cases, Nu Skin may require you to top up the price difference between your device and the new product.

Your market may offer extended warranty coverage on your locally purchased device. Please visit your local Nu Skin website or contact your local Nu Skin Support Services for more information

PATENTS

US and international patents pending.

TECHNICAL AND REGULATORY INFORMATION

ageLOC® WellSpa iO™

Model: BS2R

Maximum Output:

22V, 640 μ A, 0.414Hz

Bluetooth frequency: 2.4GHz



IPX5



IPX7






ageLOC WellSpa iO Wireless Charger

Input: 5V === min 500mA



IPX4

DEFINITION OF SYMBOLS

Symbol	Definition
	Warning
	Dispose of Batteries Properly for Your Area
	Two Levels of Shock Protection
	Charge only with specified charger
	Direct Current



©2023 NSE PRODUCTS, INC.

75 WEST CENTER STREET, PROVO, UT. 84601

DISTRIBUTED BY NU SKIN ENTERPRISES SINGAPORE PTE LTD

331 NORTH BRIDGE ROAD, #18-01 ODEON TOWERS, SINGAPORE 188720

NUSKIN.COM

36 002589 94515636/1